

MAILLEFER

Conform Fit™ Gutta-Percha

The right fit and flow you can feel





Next-generation gutta-percha. Made possible by Conform Fit[™] technology.

Gutta-percha master cones should ideally fit the shapes created by your files. Yet, traditional hand-rolled gutta-percha points have made it hard to achieve a precise match. The next level of precision performance is here — gutta-percha master cones with Conform Fit technology.



The right fit you can feel. Variable tapered cones provide a snug fit that reaches the apex for the right tug back.



Improved flow to the end. Conform Fit gutta-percha has a better flow and transfers heat up to 5 mm beyond the heat source.



Enhanced patient safety. Master cones are not made with natural rubber latex.




The simplicity of a total solution. Precisionmatched shapes and sizes remove the guesswork.



Can be used at lower working temperature. Lower temperatures reduce the risk of harm to the periodontal ligament (PDL).



Convenient tab and color markings for easy identification.

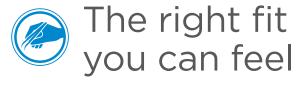



You need a consistent and confident cone fit. But finding the right tug back with hand-rolled gutta-percha points can be tough. Traditional gutta-percha may seem to fit but can actually be loose at the apex. Conform Fit™ gutta-percha cones feature a variable taper with a 40% more consistent shape than traditional gutta-percha. This creates a snug apical fit and an accurate tug back that you can feel.



Improved flow to the end

The unique gutta-percha formulation heats more consistently for better flow and allows deeper thermal penetration into the material. Heat transfers through the master cone up to 5 mm beyond the heat source—so you can be sure you have warm, flowable gutta-percha all the way to the apex. This enables the clinician to utilize a lower heat plugger temperature, thereby improving patient safety and achieving a complete fill.

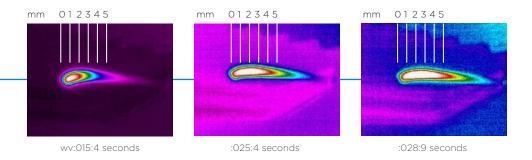

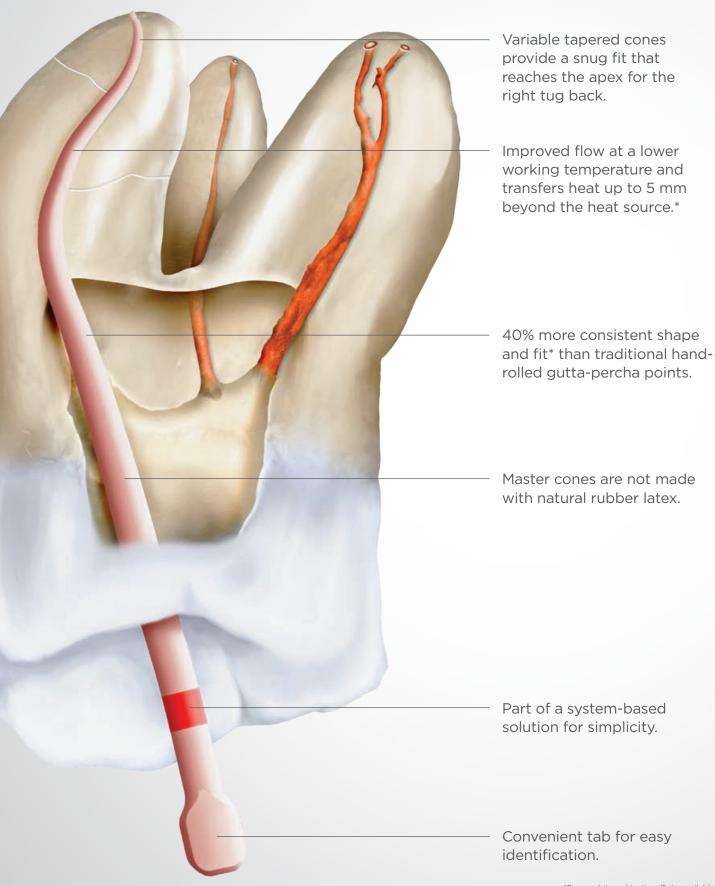



Enhanced patient safety

Conform Fit™ gutta-percha is not made with natural rubber latex. Our advanced formulation can be used at a lower working temperature, reducing the risk of harm to the periodontal ligament (PDL), which can be damaged by higher heat.

The simplicity of a total solution

Gutta-percha master cones with Conform Fit technology are part of a total system solution that takes the guesswork out of achieving a precise fit. Use them with corresponding files, gutta-percha points and paper points with confidence. The innovative tab design simplifies identification.







Conform Fit™ Gutta-Percha

Part of a total system solution

WaveOne® Gold Conform Fit™ Gutta-Percha

Color	Size	Pack	Item #
Yellow	Small	60-pack	B00WGGPF00SML
Red	Primary	60-pack	B00WGGPF00PRI
Green	Medium	60-pack	B00WGGPF00MED
■ White	Large	60-pack	B00WGGPF00LRG
Assorted	S/P/M/L	60-pack	B00WGGPF00AST





ProTaper Gold® Conform Fit™ Gutta-Percha

Size	Pack	Item #
F1	60-pack	BOOPGGPF000F1
F2	60-pack	BOOPGGPF000F2
F3	60-pack	BOOPGGPF000F3
F2/F3	60-pack	BOOPGGPF00F23
F4/F5	60-pack	BOOPGGPF00F45
	F1 F2 F3 F2/F3	F1 60-pack F2 60-pack F3 60-pack F2/F3 60-pack






ProTaper Next® Conform Fit™ Gutta-Percha

Color	Size	Pack	Item #
Red	X2	60-pack	B00PNGPF000X2
■Blue	X3	60-pack	B00PNGPF000X3
Assorted	X2/X3	60-pack	B00PNGPF00X23
Assorted	X4/X5	60-pack	BOOPNGPF00X45



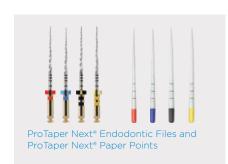

TruNatomy™ Conform Fit™ Gutta-Percha

Color	Size	Pack	Item #
Yellow	Small	60-pack	B00TNGPF00SML
Red	Prime	60-pack	B00TNGPF00PRM
■ Green	Medium	60-pack	BOOTNGPFOOMED
Assorted	S/P/M	60-pack	BOOTNGPFOOAST









